

Dios

DIOS-7300
MASSAGE CHAIR
MANUAL





7D
MASSAGE
SYSTEM





Please read and follow all safety and operating instructions listed in this user manual prior to use. Basic precautions should always be used when operating an electrical appliance.

Safety Instructions	2
Product Components	7
Chair Installation	9
Remote Controller	12
Operation Guide	14
Troubleshooting	22
Maintenance	23
Product Specifications	24

SAFETY INSTRUCTION

- Please read this manual carefully to ensure the correct installation and operation of the massage chair.
- The "Safety Instructions" section provides important safety information. Both "Warning" and "Caution" are critical safety precautions. Please pay close attention.

 Warning	This symbol indicates a warning of possible death or injury.
 Caution	This symbol indicates a caution of possible pain or damage to the body.

	This symbol indicates prohibited behavior.
	This symbol indicates instructions that must be followed for proper operation.
	This symbol indicates that self-disassembly is prohibited.
	This symbol indicates the correct operation of the cable.


- The following information should be kept for future reference.





Warning

	<ul style="list-style-type: none"> • Osteoporosis patients: Consult a doctor before use. • People with heart disease, including those with embedded pacemakers or other medical electronic devices: Consult a doctor before use. • Pregnant women or individuals with health issues: Consult a doctor before use. • Patients with a high fever, perceptual disabilities, or individuals over 80 years old are prohibited from using this machine. • Children are not allowed to use the chair. • Do not operate the chair after consuming alcohol. • Do not sleep in the chair. • Do not use the chair when your body is wet. 	<p>People who are not suitable to use the chair</p>
	<ul style="list-style-type: none"> • Ensure there are no children, pets, or objects around (under, behind, or in front of) the chair when adjusting the backrest, leg rest, or reclining position. • If you experience any discomfort during operation, stop immediately and consult a physician. • Avoid strong massage on the head. • Do not massage the belly or knees. Be cautious when massaging the neck. • Children or pets are not allowed to play on the chair, especially during operation. Standing or sitting on the backrest or armrest is prohibited. • Individuals with limited mobility or difficulty expressing themselves must be supervised continuously during use. 	<p>Accident /injury</p>
	<ul style="list-style-type: none"> • Accident/Injury Prevention: Insert the plug into a correct and suitable power socket. Ensure the plug is fully inserted to reduce the risk of short circuits or fire. • Unplug the power when not in use. Stop all massage functions and turn off the power switch before unplugging. 	<p>Fire and electric shock</p>
	<ul style="list-style-type: none"> • Do not fold the massage cushions for storage. • Do not run the machine without a load for extended periods. 	<p>Accident</p>
	<ul style="list-style-type: none"> • Do not use the power cord, plug, or socket if they are damaged or have loose connections. • Use only the rated voltage to avoid fire, electric shock, or malfunction. 	<p>Fire and electric shock</p>




Warning

	<ul style="list-style-type: none"> • If the chair is dysfunctional or broken, splashed with water, or falls into water, please have it repaired by a professional repairman or send the chair to an assigned service agent. Self-modification, disassembly, or maintenance is forbidden. • If the power cord or plug is damaged, please don't replace it yourself, please contact us or an appointed dealer for replacement. • To avoid electric shock or body injury, please don't open or disassemble any part of the chair, there is not any attached part of the chair for maintenance by customers. • Check if leather or cloth material on the back and other areas are broken or not. If broken, even if small, please stop using the chair, unplug the power, and have it repaired. 	<p>Disassembly and maintenance</p>
-----------------------------------------------------------------------------------	---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------	------------------------------------


Caution


	<ul style="list-style-type: none"> • Please operate the chair when it is flat or in a leveled place. 	<p>To avoid sudden falling.</p>
	<ul style="list-style-type: none"> • Position the chair 4 inches from the wall to allow for sufficient space to recline up and down. 	<p>To make sure the back lying.</p>
	<ul style="list-style-type: none"> • The maximum suggested operating time is 20 minutes. 5 mins for each position beginning with the slight- step massage. A continuous massage in the same position may cause uncomfortable feelings. 	<p>Too long time massage is bad for your body.</p>
	<ul style="list-style-type: none"> • Before sitting on the chair to avoid injury. Please inspect everything to confirm that the massage heads are in the correct position. 	<p>To avoid accidental injury.</p>
	<ul style="list-style-type: none"> • Please make sure that the power switch is on "OFF" position before inserting the plug into the socket. 	<p>To avoid the accident when turn on the power.</p>
	<ul style="list-style-type: none"> • Turn off all controls before removing the plug from the outlet. 	<p>To avoid the accident when turn on the power.</p>
	<ul style="list-style-type: none"> • Make sure the plug is completely inserted into the socket. 	<p>To reduce short circuits and any fire hazard.</p>
	<ul style="list-style-type: none"> • Check whether there are any damages on the wire and make sure there is no water or rubbish on the wire. 	<p>To avoid electric shock in any leakage or breakdown.</p>
	<ul style="list-style-type: none"> • Don't use or store the chair in a wet environment like near the swimming pool, or bathroom. • Keep the chair and remote control dry at all times. • Don't use the chair when your body is wet. 	<p>To avoid leakage or electric shock.</p>

Caution

	<ul style="list-style-type: none"> • Don't press with your hands or step on the massage heads with your feet. • Don't put your hands between the working massage heads. • Don't put your hands or feet in the gap between mechanical parts. • Don't drop or insert anything into the slot of the massage chair. • Don't stand or sit or put anything on the backrest, armrest, or legrest. • Don't stand on the chair during the operation. Don't move the chair during the operation. • Don't pull out the plug or turn off the power during operation. • Don't put your hands or head between leg-rest or beneath. • Don't use the massage chair while wearing hair ornaments. • Don't put hard things in your pocket when taking buttocks and huckle massage. • Don't sit or press on the backrest when the chair is in a reclining position. Don't use it on bare skin or wearing thin clothes, thin clothes may cause a better massage effect, But it also may cause irritation to the skin or cryogenic burns. 	<p>To avoid accidents or injury</p>
	<ul style="list-style-type: none"> • Keep the chair away from heat, wet, sharp things and corrosive, inflammable, explosive environment. • Using an electric blanket or other heating products on the chair is forbidden. • Don't put in the space of a serious or corrosive gas environment. • Don't put the chair in small space places or use it in bad air circulation places. <hr/> <ul style="list-style-type: none"> • Keep the chair away from the stove or other heating products. • Don't put the chair directly under the sun or other high temperature places. <hr/> <ul style="list-style-type: none"> • Don't use the chair outdoors. 	<p>To avoid fade, discolor on upholstery or hardened, deteriorated on PU.</p>
	<ul style="list-style-type: none"> • Make sure to install the wire in the correct way(In-correct installation may cause electric shock. If you don't know how to install, please find a professional electrician) • Don't use a plug adapter (If the enclosed plug doesn't match with socket, please don't change the plug by yourself and find a professional electrician to install a matched socket.) 	<p>To avoid electric shock.</p>


 **Caution**

	<ul style="list-style-type: none"> • Unplug the power before cleaning the chair. • Unplug the power after using the chair. 	To avoid electric shock or fire hazards.
	<ul style="list-style-type: none"> • Unplug the power when any damage on the power. 	
	<ul style="list-style-type: none"> • Hold the plug to unplug the power but to not try to pull the wire to unplug the power. 	To avoid accidents
	<ul style="list-style-type: none"> • Please turn off the power switch and unplug the power after massage. 	
	<ul style="list-style-type: none"> • Unplug the power immediately when electricity goes off. 	


-  **Warning :**
- 1.Do not use it when gets wet.
 - 2.The disabled people who need supervision.
 - 3.The chair has heat function, people who are insensitive to heat must pay high attention.

| Grounding instruction |

The chair is equipped with a cord having a grounding conductor. The plug must be plugged into an appropriate socket that is properly installed and matched with the plug.

 <p>The earth wire is not allowed to be in these places</p>	<p>Gas pipe: probably cause explosions or fire. Phone wire and lightning rod: probably cause electric shock and fire when thundering. Water pipe: earth wire is useless where there are plastics</p>
------------------------------------------------------------------------------------------------------------------------------------------------	--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------


| Main body diagram (front) |

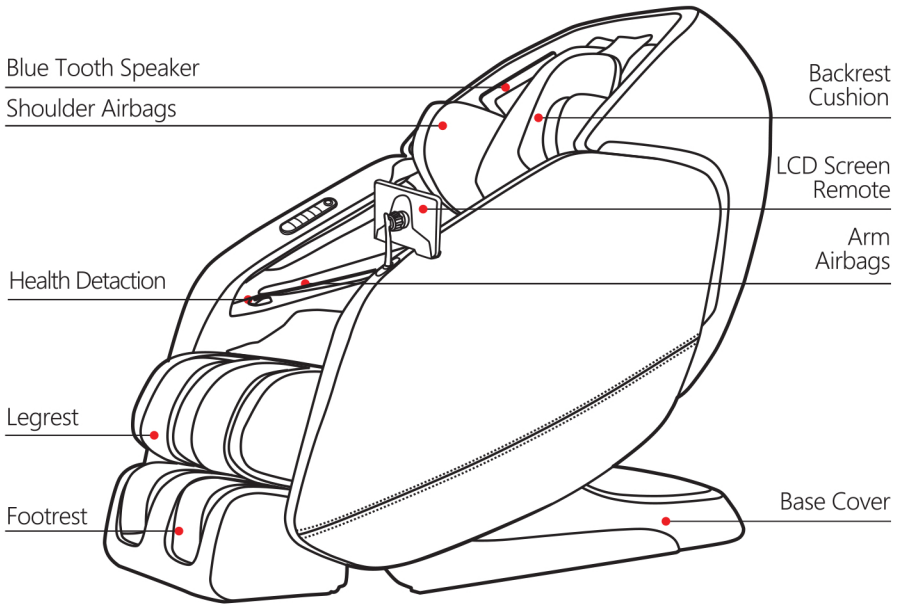


Built in Premium Bluetooth Speakers



| Main body diagram (rear) |

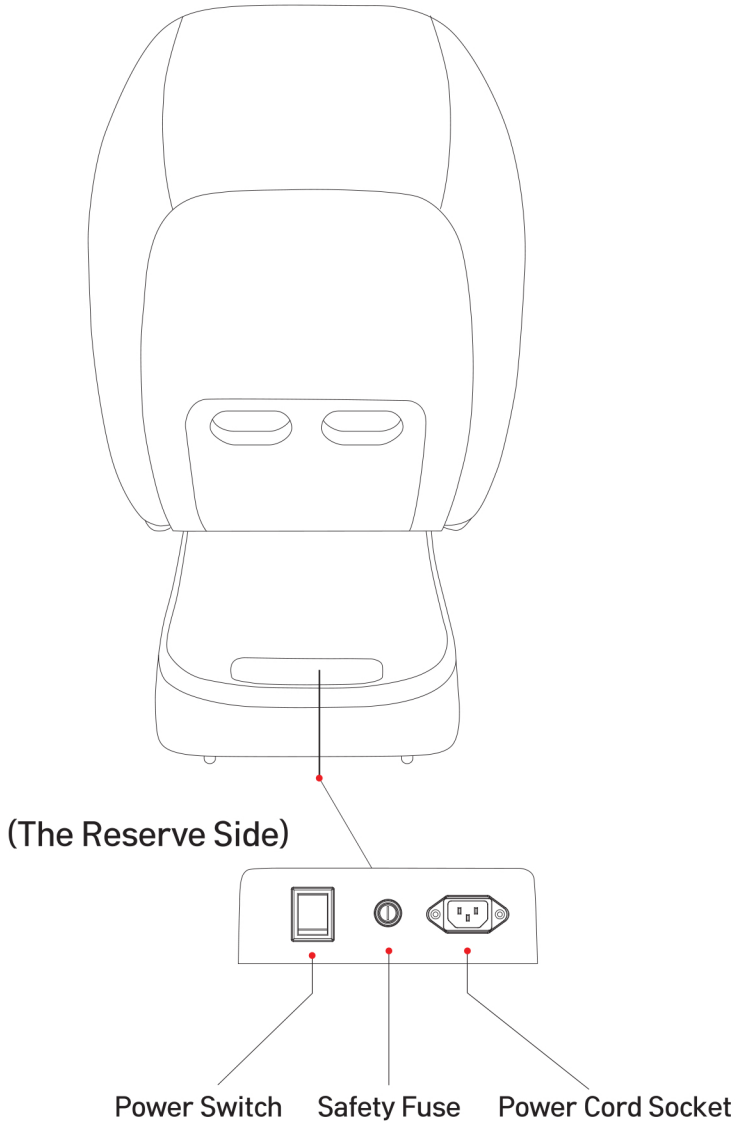


1. Install leg-rest

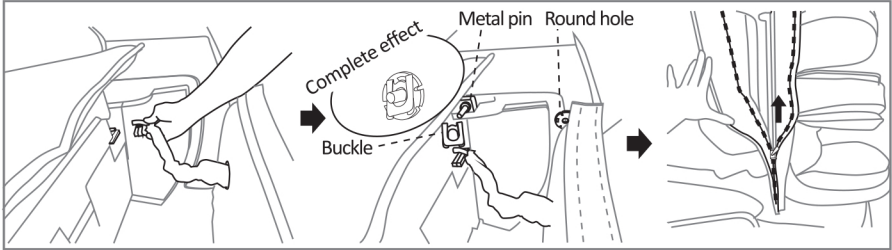

1. First connect the trachea of the leg with the bus on the massage chair.
2. Lift the leg and insert the round hole in the upper right corner of the leg into the metal pin sha of the massage chair, and then insert it with a buckle to fix it.
3. When installing the le side, first move the whole leg to the right, align the round hole in the upper le corner and the pin of the massage chair, move the lower leg to the left, and then insert the buckle to fix it.
4. Pull the holster on the leg to the zipper of the massage chair.
5. Finish the leather case.

2. Install the remote controller

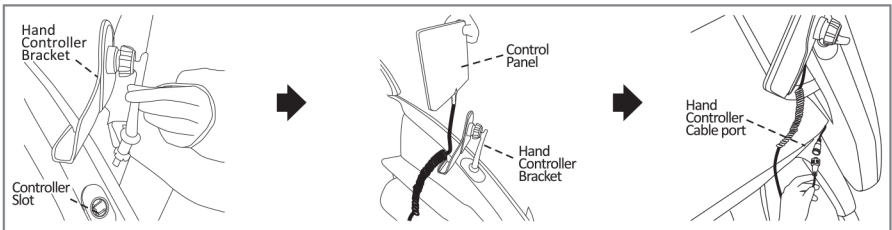


1. Firstly, snap the hand controller bracket into the hand controller slot according to the direc on of the slot.
2. Place the hand controller on the controller bracket.
Complete effect
3. Lastly insert handle controller cable port is docked with the port of the massage chair seat.

CHAIR INSTALLATION


| 3. Before use |

Use the controller to place the backrest frame to the extreme position, raise the calf to the highest position and raise it to the longest position. Make sure there is sufficient space around the massage chair before the operation.

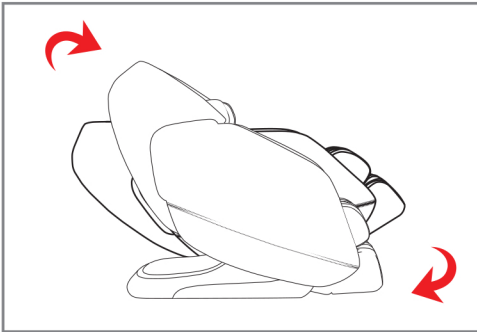

Use the controller to place the backrest frame to the extreme position, raise the calf to the highest position and raise it to the longest position. Make sure there is sufficient space around the massage chair before the operation.


4. Place the massage chair in an open area |

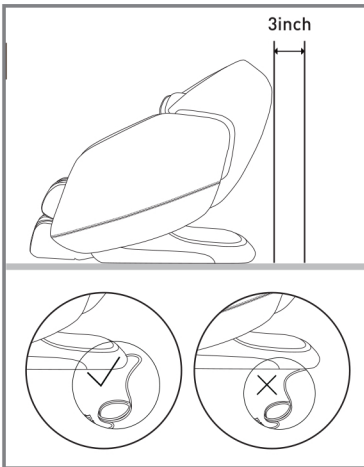

- Allow sufficient space around the chair for recline purposes.
- 3 inch distance for the back side and 19.7 inches in front of the footrest.
- * Do not put the power cord under the relaxed chair!

| 5. Moving the massage chair |

After completing the installation, put the chair in an incline / upright position to move the chair easily.

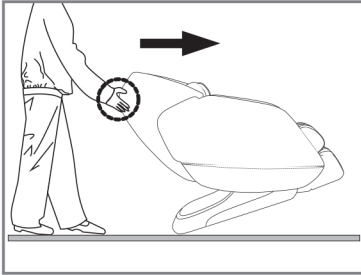


- As illustrated, tilt the chair back onto its wheels at 45° angle. Slowly push and maneuver the massage chair to the desired location. Lower the massage chair slowly and smoothly until it is fully flat on the ground.

CAUTION

- * Do not move the massage chair when occupied or in use.
- * Moving the massage chair may cause damage to the floors. Floor protection is recommended (e.g. rugs or carpets).
- * Do not drop or suddenly let go of the chair when moving as it may damage the chair.
- * Two persons are recommended when moving and extra care when handling the massage chair.
- * Make sure the power plug & remote is disconnected to prevent damage. The warranty will not cover the power cord & remote if the chair is moved with the plugs connected.

| Quick Control |

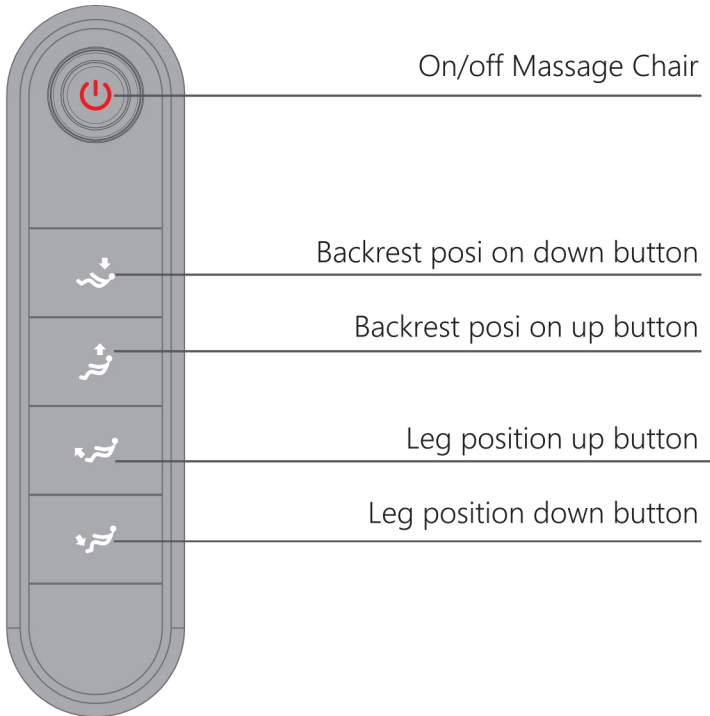

| Remote Controller |




- Power switch button : Turn the whole machine massage function on/off.
- In the powered on state, press this button to turn on the massage chair, then select other buttons (such as automatic programs, massage techniques, etc.) for a massage experience.

OPERATION GUIDE

1.Remote contoller interface and description

- Open the remote controller and the main menu will be displayed. Please select the program you want.

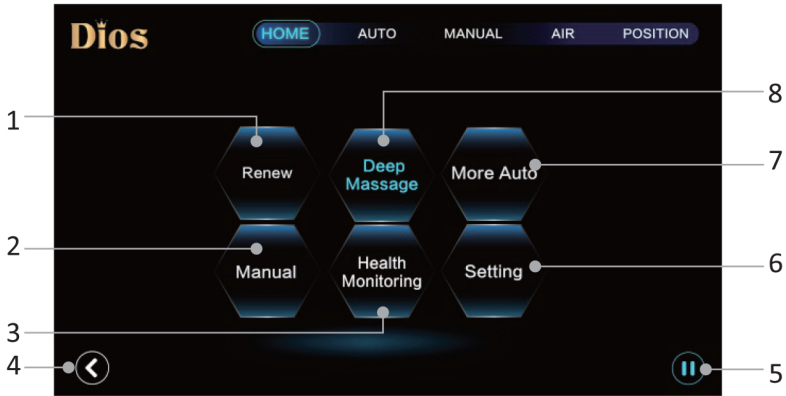

1.Renew

2.Manual modes

3.Health Monitoring

4.Return


5.Pause

6.Settings

7.More automatic modes

8.Deep massage

2.Settings

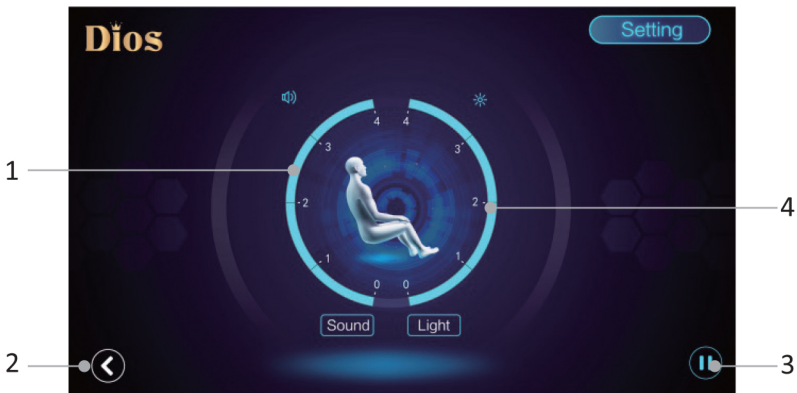

1.Sound adjustment area

2.Return


3.Pause

4.Screen brightness adjustment area

3. Auto modes




- There are 24 kinds of all-round unique automatic massage modes in the corresponding massage area.

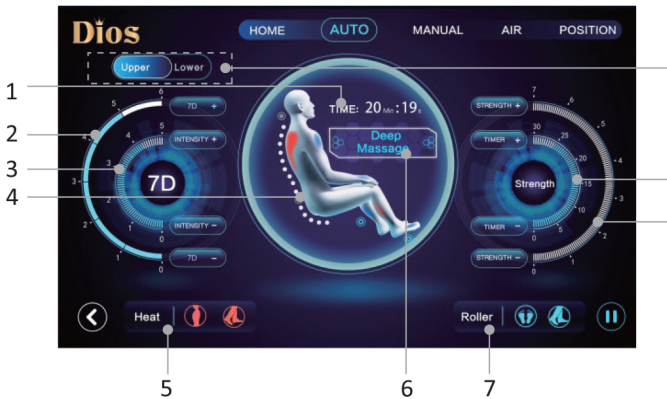

1. Display the remaining massage me
2. 7D massage intensity adjustment area
3. Airbags intensity adjustment zone
4. Movement/airbag/hea ng/roller massage area display
5. Hea ng site selec on
6. Automa c program name display
7. Sole roller loca on selec on
8. Movement speed/strength adjustment
9. Massage me adjustment
10. Manual massage of the upper or lower body


Statement: Dual movement opera on rules follow the lower movement to avoid the upper movement opera on when the distance between the two movements triggers an-collision safety logic, the lower movement will give priority to traveling down to avoid a crash with the upper movement, so it will be a li le different from the set mode of opera on, it is a normal phenomenon.

4. Detection

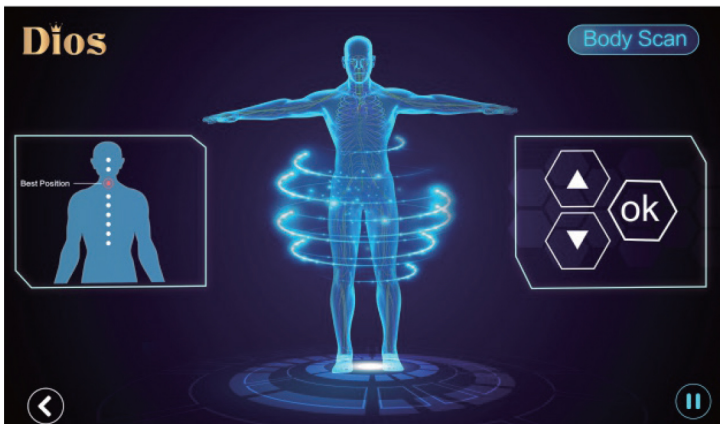


1. Press any of the automatic modes, the latest smart body detection will be turned on (only automatically detected once after booting).
2. After the inspection is completed, within 8 seconds, the user can choose to fine-tune the shoulder position according to his needs. You can also press OK directly to confirm the shoulder position.

5. Manual modes

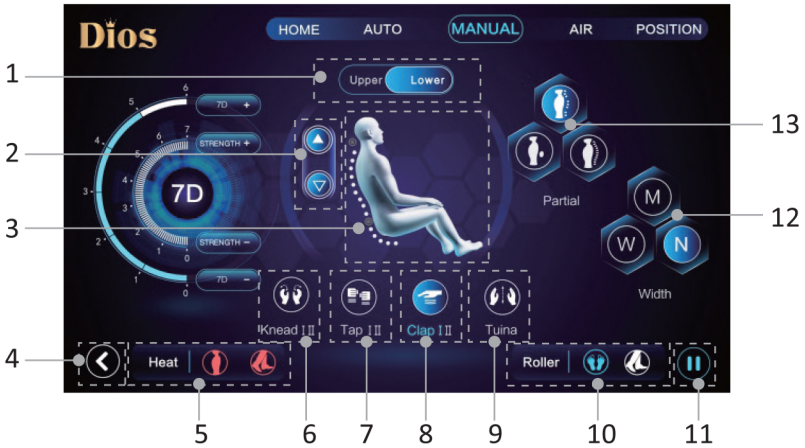


1. Manual massage of the upper or lower body
2. Adjustment of the up and down position of the massage movement
3. Massage movement area display
4. Return
5. Heating site selection
6. Kneading massage(2 levels adjustable)
7. Tapping massage(2 levels adjustable)
8. Clapping massage(2 levels adjustable)
9. Tuina Massage
10. Selection area of sole roller part
11. Pause
12. Massage width adjustment area(not effective when kneading is on)
13. Massage interval adjustment area


6. Airbags massage

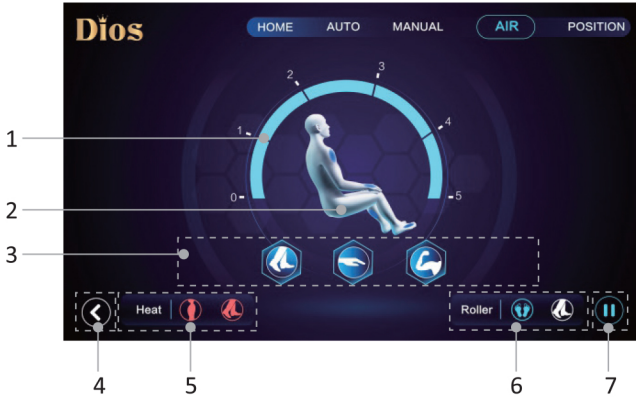

- | | |
|---------------------------------------------|---------------------------------------|
| 1. Airbag massage intensity adjustment area | 2. Airbag massage area display |
| 3. Airbag massage area selection area | 4. Return |
| 5. Heating area selection area | 6. Selection area of sole roller part |
| 7. Pause | |

7. Health monitoring

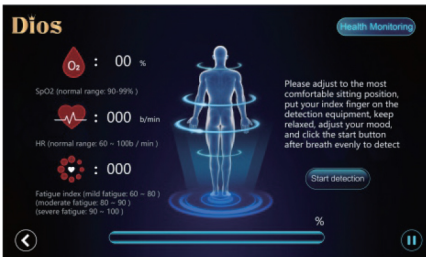


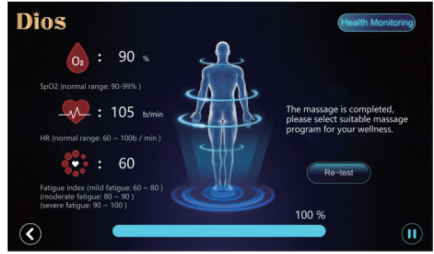

Put your index finger on the detect on device. After the detection is completed, your blood oxygen level, heart rate, and fatigue index will be displayed.

- Note: This data is for reference only and should not be used for pathological assessment.

8. Sitting adjustment

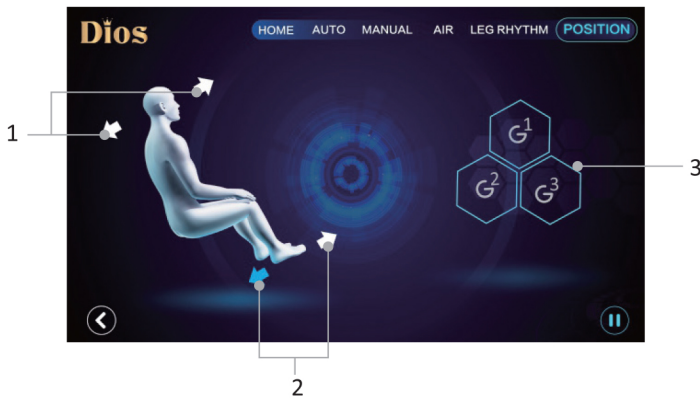

1. Backrest down and up
2. Footrest raised and lowered
3. Zero-gravity position area: three levels option


9. Pairing massage chair with bluetooth music

- a. Pairing: turn on bluetooth device in the cellphone or ipad.
Search for Dios 7300.
- b. Playing music: after pairing use the control device in the cellphone or ipad to manage.

Note: disconnect bluetooth when the massage chair is not in use




| 3 After use |

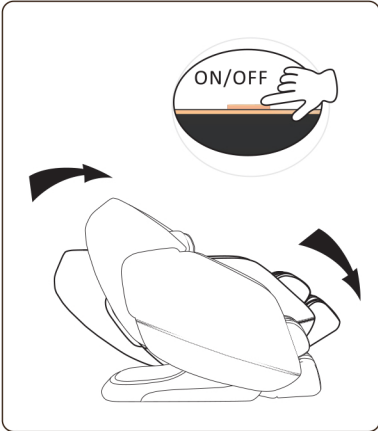

1. After the massage program is completed return the massage chair to an upright position by turning off the massage chair.

Please check the surroundings prior to adjusting the recline angle, paying special attention to children or pets that may be in close proximity.

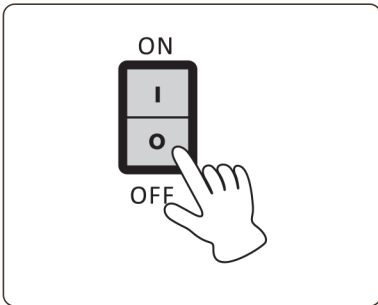

2. Press the power button to power off the massage chair and return the remote controller to the storage area in the arm pocket.

3. Turn off the main power switch when you are not using the chair for a long period. Remove Plug from the power outlet to avoid accidents by kids and pets.

4. Do not operate with wet hands.

TROUBLE SHOOTING

If the sounds or movements are not normal, please turn off the massage chair and remove the plug from the electric outlet. Please contact the customer service department for further assistance.


1.855.646.9500 / 1.714.823.3338, service@kahunachair.com

Problem	Solution
The massage chair is connected to power but is not functioning.	Check to see if the power plug is inserted or if the power outlet is functioning properly.
The remote controller is responding, but the massage chair is not functioning.	The massage chair may have automatically shut down due to overheating protection devices. Please turn off the power and allow the chair to cool off for 30 to 50 minutes.
Both the remote controller and massage chair are not responding.	The massage chair automatically turns into "full stop" mode when a massage program is completed. Please press the power button again to turn off the chair.




KAHUNA CHAIR PRODUCT WARRANTY REGISTRATION
WWW.KAHUNASERVICE.COM

| Use a dry cloth to clean |

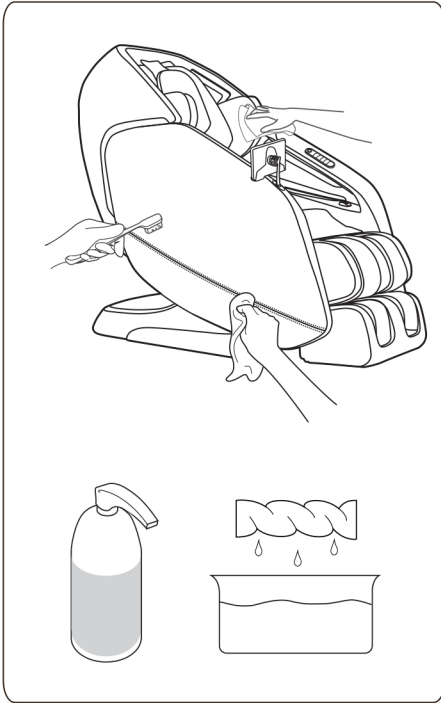

Turn off and pull out the plug before cleaning the massage chair, do not touch the plug with wet hands.

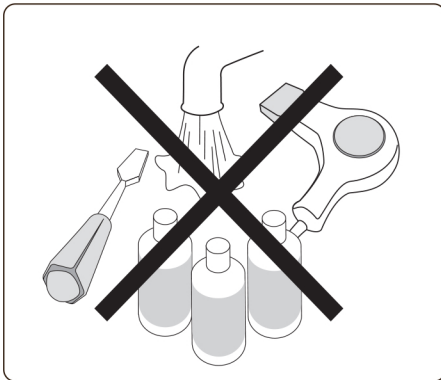

Prohibit using liquid cleaning products like water, alcohol, or other cleaning products to maintain a good quality.

PRODUCT SPECIFICATIONS

Description	Specifications
Product Name	Dios Massage Chair
Model No.	Dios-7300
Rated Voltage	AC 110V~220V
Rated Frequency	60Hz
Rated Power Consumption	200 W
Dimension (LxWxH) Inches	58.3 x 32.7 x 50.4 inches
Weight (LBS)	Net Weight : 229 lbs Gross Weight : 242.5 lbs



Product Massage Chair

Model DIOS-7300


Developed by RJ Distribution INC

11125 S Knott Ave, Suite D

Cypress, CA 90630

Tel 1-855-646-9500

Email INFO@KAHUNACHAIR.COM

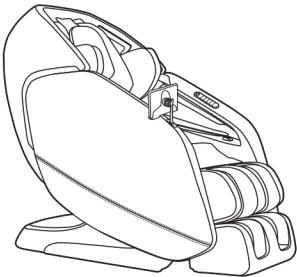

You must register warranty at www.kahunaservice.com

INSTRUCTION MANUAL
MODEL:DIOS-7300

Please read this user manual thoroughly prior to operating the massage chair. Store this booklet for future reference.